



# Filtek™

Bulk Fill  
Flowable Restorative

## The bulk fill flowable that's easier to inject.

3M™ Filtek™ Bulk Fill Flowable Restorative is available in a new ergonomic-designed syringe for improved comfort and easier injecting. This helps give you more control and visibility for precise placement in deep restorations.

### Helps reduce bubbles and waste.

Virtually no bubbles or material "run-on" during dispensing gives you better control.

### No more guessing.

The green syringe barrel indicates the remaining material volume. When your syringe is empty, only the white plunger shows.


### Easy to hold and inject.

Triangle finger plate and plunger help to make it quick and easy to personalize your grip.



### Bendable cannula.

Easy-to-bend without kinking for better access.

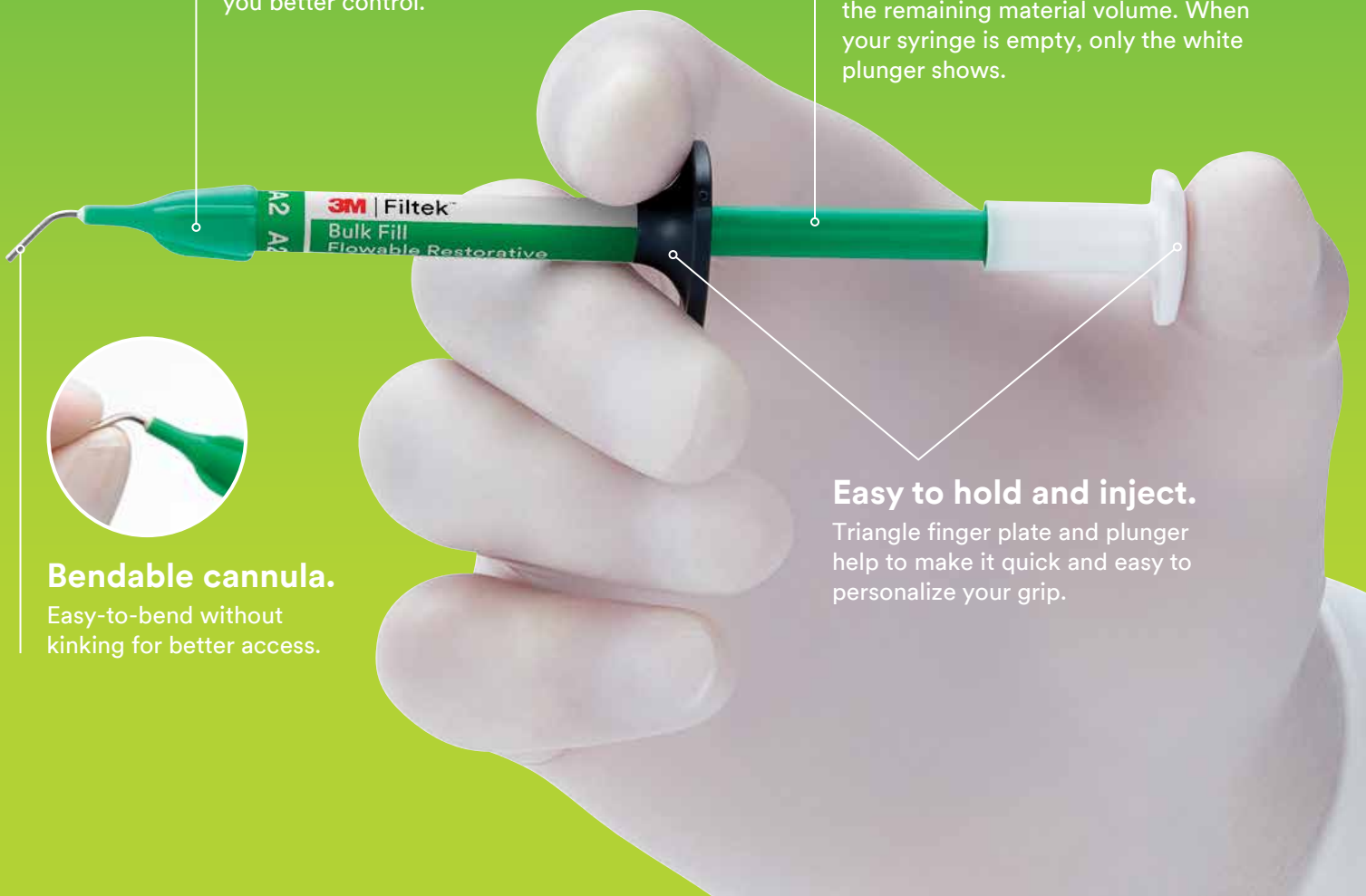

## Key benefits

- Virtually no bubbles or material “run-on”
- Easy-to-bend cannula for improved access
- Ergonomic-design is easy to hold and inject
- Injectable, one-step placement up to 4 mm
- Excellent adaptation and low shrinkage stress

## Can be warmed for injection moulding.

No higher than 70°C (158°F), no longer than 1 hour; up to 25 times.



## Indications

(Incremental up to 4 mm)



Base  
(Class I – II)



Liner  
(Class I – II)



Minimally Invasive  
Dentistry (MID)<sup>1</sup>



Class III



Class V



Repair of Indirect  
(small defects)



Repair<sup>2</sup>



Undercut  
Block-out



Pit & Fissure  
Sealant



Core  
Build-up<sup>3</sup>

Refer to Instructions  
for Use for full list of  
indications.

## Posterior restoration made simple.



Prepare the tooth and apply an adhesive e.g. 3M™ Scotchbond™ Universal Adhesive.



Dispense up to 4 mm of 3M™ Filtek™ Bulk Fill Flowable Restorative into the prepared cavity and cure, e.g. with 3M™ Elipar™ DeepCure-S LED Curing Light.



Place a 2 mm (minimum) capping layer of an aesthetic composite, such as 3M™ Filtek™ Universal Restorative, and cure.

Clinical photos courtesy of Dr. Reich, Seefeld, Germany.

## Ordering Information



### 3M™ Filtek™ Bulk Fill Flowable Restorative Refills

Syringe refills include: 2 – 2 g syringes; 20 dispensing tips.

Shade	Item No.
A1	4862A1
A2	4862A2
A3	4862A3
Universal	4862U

Item No.	Product Info.
3740T	3M™ Filtek™ Bulk Fill Flowable Restorative Dispensing Tips, 19-gauge, Green (20)



Available in Canada from

**3M Canada**  
**3M Oral Care**  
300 Tartan Drive  
London, ON N5V 4M9  
Canada  
1-888-363-3685  
3M.ca/dental

**3M ESPE Dental Products**  
2510 Conway Avenue  
St. Paul, MN 55144  
USA  
1-800-634-2249

<sup>1</sup> Including small, non stress-bearing occlusal restorations.

<sup>2</sup> Repair of acrylic temporary and composite resin.

<sup>3</sup> Where at least half the coronal tooth structure is remaining to provide structural support for the crown.